



REF 10215
β2 GP1 IgGAM

Instruction Manual

Table of Contents

1	Intended Use	1
2	Clinical Application and Principle of the Assay	1
3	Kit Contents	2
4	Storage and Shelf Life	2
5	Precautions of Use.....	3
6	Sample Collection, Handling and Storage.....	4
7	Assay Procedure	5
8	Quantitative Interpretation.....	8
9	Technical Data.....	9
10	Performance Data.....	9
11	Literature	10



AIDA GmbH
Dr.-Karl-Aschoff-Straße 9
55543 Bad Kreuznach
Germany
Phone: +49 671 92065090
Fax: +49 671 92065091
Website: www.aida-diagnostics.com
Mail: info@aida-diagnostics.com

 autoimmune diagnostic assays	Product Ref.	10215
	Product Desc.	β2 GP1 IgGAM
	Manual Rev. No.	006: 2024-02-15

1 Intended Use

β2 GP1 IgGAM is a solid phase enzyme immunoassay employing native β2-glycoprotein I highly purified from human plasma for the combined quantitative detection of IgA, IgG and IgM antibodies against β2-glycoprotein I in human serum. Anti-β2-glycoprotein I antibodies recognize specific epitopes on human β2-glycoprotein I which are expressed only when β2-glycoprotein I interacts with lipid membranes or when absorbed to other surfaces (e.g. microtiter plate). The assay is an aid in the diagnosis and risk of primary and secondary antiphospholipid syndrome (APS).

2 Clinical Application and Principle of the Assay

Antibodies against β2-glycoprotein I belong to the group of anti-phospholipid antibodies mainly targeted against complexes composed of negatively charged phospholipids (cardiolipin e.g) and plasma proteins like β2- glycoprotein I, prothrombin, protein C or protein S.

Reactivity against isolated β2- glycoprotein I is found, too. Thus β2-glycoprotein I is discussed to be an autoantigen on its own. β2-glycoprotein I, also called apolipoprotein H, is a 50 kDa β2 globuline which is associated in vivo with lipoprotein, platelets and phospholipids and which seems to inhibit the intrinsic coagulation pathway, the prothrombinase activity and the ADP-dependent platelet aggregation. Anti-phospholipid antibodies are frequently found in sera of patients with systemic lupus erythematosus and related diseases and are typical for the secondary development of an antiphospholipid syndrome (APS). Whilst, anti-phospholipid antibodies in patients with no other autoimmune diseases characterize a primary APS.

Many studies have shown a correlation between these autoantibodies and an enhanced incidence of thrombosis, thrombocytopenia and habitual abortions (as a consequence of placental infarction). The exact mechanism how pathogenic anti-phospholipid antibodies induce thrombosis has not been revealed yet.

Principle of the test

Serum samples diluted 1:101 are incubated in the microplates coated with the specific antigen. Patient's antibodies, if present in the specimen, bind to the antigen. The unbound fraction is washed off in the following step. Afterwards anti-human immunoglobulins conjugated to horseradish peroxidase (conjugate) are incubated and react with the antigen-antibody complex of the samples in the microplates. Unbound conjugate is washed off in the following step. Addition of TMB-substrate generates an enzymatic colorimetric (blue) reaction, which is stopped by diluted acid (color changes to yellow). The intensity of color formation from the chromogen is a function of the amount of conjugate bound to the antigen-antibody complex and this is proportional to the initial concentration of the respective antibodies in the patient sample.

	Product Ref.	10215
	Product Desc.	β2 GP1 IgGAM
	Manual Rev. No.	006: 2024-02-15

3 Kit Contents

TO BE RECONSTITUTED				
Item	Quantity	Cap color	Solution color	Description / Contents
Sample Buffer (5x)	1 x 20ml	White	Yellow	5 x concentrated Tris, sodium chloride (NaCl), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Wash Buffer (50x)	1 x 20ml	White	Green	50 x concentrated Tris, NaCl, Tween 20, sodium azide < 0.1% (preservative)
READY TO USE				
Item	Quantity	Cap color	Solution color	Description / Contents
Negative Control	1 x 1.5ml	Green	Colorless	Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Positive Control	1 x 1.5ml	Red	Yellow	Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Calibrators	6 x 1.5ml	White	Yellow *	Concentration of each calibrator: 0, 3, 10, 30, 100, 300 U/ml. Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Conjugate, IgA/G/M	1 x 15ml	White	Red	Containing: Anti-human immunoglobulins conjugated to horseradish peroxidase, bovine serum albumin (BSA)
TMB Substrate	1 x 15ml	Black	Colorless	Stabilized tetramethylbenzidine and hydrogen peroxide (TMB/H ₂ O ₂)
Stop Solution	1 x 15ml	White	Colorless	1M Hydrochloric Acid
Microtiter plate	12 x 8 well strips	N/A	N/A	With breakaway microwells. Refer to paragraph 1 for coating.
* Color increasing with concentration				
MATERIALS REQUIRED, BUT NOT PROVIDED				
Microtiter plate reader 450 nm reading filter and recommended 620 nm reference filter (600-690 nm). Glass ware (cylinder 100-1000ml), test tubes for dilutions. Vortex mixer, precision pipettes (10, 100, 200, 500, 1000 µl) or adjustable multipipette (100-1000µl). Microplate washing device (300 µl repeating or multichannel pipette or automated system), adsorbent paper. Our tests are designed to be used with purified water according to the definition of the United States Pharmacopeia (USP 26 - NF 21) and the European Pharmacopeia (Eur.Ph. 4th ed.).				

4 Storage and Shelf Life

Store all reagents and the microplate at 2-8°C/35.6-46.4°F, in their original containers. Once prepared, reconstituted solutions are stable at 2-8°C/35.6-46.4°F for at least 1 month. Reagents and the microplate shall be used within the expiry date indicated on each component, only. Avoid intense exposure of TMB solution to light. Store microplates in designated foil, including the desiccant, and seal tightly.

 autoimmune diagnostic assays	Product Ref.	10215
	Product Desc.	β2 GP1 IgGAM
	Manual Rev. No.	006: 2024-02-15

5 Precautions of Use

5.1 Health hazard data

THIS PRODUCT IS FOR IN VITRO DIAGNOSTIC USE ONLY. Thus, only staff trained and specially advised in methods of in vitro diagnostics may perform the kit. Although this product is not considered particularly toxic or dangerous in conditions of the intended use, refer to the following for maximum safety:

Recommendations and precautions

This kit contains potentially hazardous components. Though kit reagents are not classified being irritant to eyes and skin we recommend to avoid contact with eyes and skin and wear disposable gloves.

WARNING! Calibrators, Controls and Buffers contain sodium azide (NaN_3) as a preservative. NaN_3 may be toxic if ingested or adsorbed by skin or eyes. NaN_3 may react with lead and copper plumbing to form highly explosive metal azides. On disposal, flush with a large volume of water to prevent azide build-up. Please refer to decontamination procedures as outlined by CDC or other local/national guidelines.

Do not smoke, eat or drink when manipulating the kit. Do not pipette by mouth.

All human source material used for some reagents of this kit (controls, standards e.g.) has been tested by approved methods and found negative for HbsAg, Hepatitis C and HIV 1. However, no test can guarantee the absence of viral agents in such material completely. Thus handle kit controls, standards and patient samples as if capable of transmitting infectious diseases and according to national requirements.

The kit contains material of animal origin as stated in the table of contents, handle according to national requirements.

 autoimmune diagnostic assays	Product Ref.	10215
	Product Desc.	β2 GP1 IgGAM
	Manual Rev. No.	006: 2024-02-15

5.2 General directions for use

In case that the product information, including the labeling, is defective or incorrect please contact the manufacturer or the supplier of the test kit.

Do not mix or substitute Controls, Calibrators, Conjugates or microplates from different lot numbers. This may lead to variations in the results.

Allow all components to reach room temperature (20-32°C/68-89.6°F) before use, mix well and follow the recommended incubation scheme for an optimum performance of the test.

Incubation: We recommend test performance at 30°C/86°F for automated systems.

Never expose components to higher temperature than 37°C/ 98.6°F.

Always pipette substrate solution with brand new tips only. Protect this reagent from light. Never pipette conjugate with tips used with other reagents prior.

A definite clinical diagnosis should not be based on the results of the performed test only, but should be made by the physician after all clinical and laboratory findings have been evaluated. The diagnosis is to be verified using different diagnostic methods.

6 Sample Collection, Handling and Storage

Use preferentially freshly collected serum samples. Blood withdrawal must follow national requirements. Do not use icteric, lipemic, hemolysed or bacterially contaminated samples. Sera with particles should be cleared by low speed centrifugation (<1000 x g). Blood samples should be collected in clean, dry and empty tubes.

After separation, the serum samples should be used during the first 8h, respectively stored tightly closed at 2-8°C/35.6-46.4°F up to 48h, or frozen at -20°C/-4°F for longer periods.

 autoimmune diagnostic assays	Product Ref.	10215
	Product Desc.	β2 GP1 IgGAM
	Manual Rev. No.	006: 2024-02-15

7 Assay Procedure

7.1 Preparations prior to starting

Dilute concentrated reagents:

Dilute the concentrated sample buffer 1:5 with distilled water (e.g. 20 ml plus 80 ml).

Dilute the concentrated wash buffer 1:50 with distilled water (e.g. 20 ml plus 980 ml).

To avoid mistakes we suggest to mark the cap of the different calibrators.

Samples:

Dilute serum samples 1:101 with sample buffer (1x)

e.g. 1000 µl sample buffer (1x) + 10 µl serum. Mix well !

Washing:

Prepare 20 ml of diluted wash buffer (1x) per 8 wells or 200 ml for 96 wells

e.g. 4 ml concentrate plus 196 ml distilled water.

Automated washing:

Consider excess volumes required for setting up the instrument and dead volume of robot pipette.

Manual washing:

Discard liquid from wells by inverting the plate. Knock the microwell frame with wells downside vigorously on clean adsorbent paper. Pipette 300 µl of diluted wash buffer into each well, wait for 20 seconds. Repeat the whole procedure twice again.

Microplates:

Calculate the number of wells required for the test. Remove unused wells from the frame, replace and store in the provided plastic bag, together with desiccant, seal tightly (2-8°C/35.6-46.4°F).

7.2 Pipetting Scheme

We suggest pipetting calibrators, controls and samples as follows:

For *QUANTITATIVE* interpretation

	1	2	3	4...
A	Cal A	Cal E	P1	
B	Cal A	Cal E	P1	
C	Cal B	Cal F	P2	
D	Cal B	Cal F	P2	
E	Cal C	PC	P3	
F	Cal C	PC	P3	
G	Cal D	NC	...	
H	Cal D	NC	...	

CalA: calibrator A

CalB: calibrator B

CalC: calibrator C

CalD: calibrator D

CalE: calibrator E

CalF: calibrator F

PC: positive control



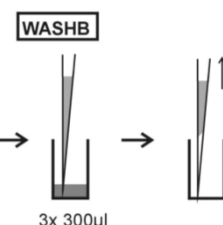
NC: negative control


P1: patient 1


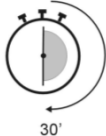
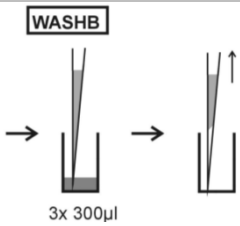
P2: patient 2


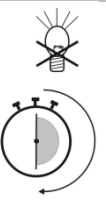
P3: patient 3

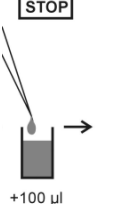

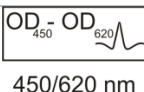
7.3 Test Steps

Step	Description
1.	Ensure preparations from step 7.1 above have been carried out prior to pipetting.
2.	Use the following steps in accordance with quantitative interpretation results desired:
CONTROLS & SAMPLES	
3.	 <p>Pipette into the designated wells as described in chapter 7.2 above, 100 µl of either: Calibrators (CAL.A to CAL.F) and 100 µl of each of the following:</p> <ul style="list-style-type: none"> Negative control (NC) and Positive control (PC), and Patients diluted serum (P1, P2...)
4.	 <p>Incubate for 30 minutes at 20-32°C/68-89.6°F.</p>
5.	 <p>Wash 3x with 300 µl washing buffer (diluted 1:50).</p>

<div> autoimmune diagnostic assays</div>		Product Ref.	10215
		Product Desc.	β2 GP1 IgGAM
		Manual Rev. No.	006: 2024-02-15

CONJUGATE		
6.	<div><div>CONJ</div></div>	Pipette 100 µl conjugate into each well.
7.	<div></div>	Incubate for 30 minutes at 20-32°C/68-89.6°F.
8.	<div><div>WASHB</div></div>	Wash 3x with 300 µl washing buffer (diluted 1:50).

SUBSTRATE		
9.	<div><div>SUB</div></div>	Pipette 100 µl TMB substrate into each well.
10.	<div></div>	Incubate for 30 minutes at 20-32°C/68-89.6°F, protected from intense light.

STOP		
11.	<div><div>STOP</div></div>	Pipette 100 µl stop solution into each well, using the same order as pipetting the substrate.
12.	<div></div>	Incubate 5 minutes minimum.
13.		Agitate plate carefully for 5 sec.
14.	<div><div>OD₄₅₀ - OD₆₂₀</div><div>450/620 nm</div></div>	Read absorbance at 450 nm (recommended 450/620 nm) within 30 minutes.

8 Quantitative Interpretation

For **quantitative interpretation** establish the standard curve by plotting the **optical density (OD) of each calibrator (y-axis)** with respect to the corresponding concentration values in U/ml (x-axis). For best results we recommend log/lin coordinates and 4-Parameter Fit. From the OD of each sample, read the corresponding antibody concentrations expressed in U/ml.

Normal Range	Equivocal Range	Positive Results
< 16 U/ml	16 - 24 U/ml	>24 U/ml

Example of a standard curve

Do NOT use this example for interpreting patient's result

Calibrators IgG/A/M	OD 450/620 nm	CV % (Variation)
0 U/ml	0.041	1.7
3 U/ml	0.132	0.0
10 U/ml	0.280	2.6
30 U/ml	0.584	2.1
100 U/ml	1.211	0.0
300 U/ml	2.042	0.6

Example of calculation

Patient	Replicate (OD)	Mean (OD)	Result (U/ml)
P 01	0.756/0.757	0.757	45.9
P 02	1.344/1.352	1.348	123.3

Samples above the highest calibrator range should be reported as >Max. They should be diluted as appropriate and re-assayed. Samples below calibrator range should be reported as < Min.

For lot specific data, see enclosed quality control leaflet. Medical laboratories might perform an in-house quality control by using own controls and/or internal pooled sera, as foreseen by national regulations.

Each laboratory should establish its own normal range based upon its own techniques, controls, equipment and patient population according to their own established procedures.

In case that the values of the controls do not meet the criteria the test is invalid and has to be repeated.

The following technical issues should be verified: Expiration dates of (prepared) reagents, storage conditions, pipettes, devices, photometer, incubation conditions and washing methods.

If the items tested show aberrant values or any kind of deviation or that the validation criteria are not met without explicable cause please contact the manufacturer or the supplier of the test kit.

 autoimmune diagnostic assays	Product Ref.	10215
	Product Desc.	β2 GP1 IgGAM
	Manual Rev. No.	006: 2024-02-15

9 Technical Data

Sample material:	serum
Sample volume:	10 µl of sample diluted 1:101 with 1x sample buffer
Total incubation time:	90 minutes at 20-32°C/68-89.6°F
Calibration range:	0-300 U/ml
Analytical sensitivity:	1.0 U/ml
Storage:	at 2-8°C/35.6-46.4°F use original vials only.
Number of determinations:	96 tests

10 Performance Data

10.1 Analytical sensitivity

Testing sample buffer 30 times on **β2 GP1 IgGAM** gave an analytical sensitivity of 1.0 U/ml.

10.2 Specificity and sensitivity

The microplate is coated with highly purified native human β2-Glycoprotein I. No crossreactivities to other autoantigens have been found. The β2 GP1 test exhibits a diagnostic specificity of 100%. The β2 GP1 test exhibits a diagnostic sensitivity of 47%.

10.3 Linearity

Chosen sera have been tested with this kit and found to dilute linearly. However, due to the heterogeneous nature of human autoantibodies there might be samples that do not follow this rule.

Sample No.	Dilution Factor	Measured (U/ml)	Expected (U/ml)	Recovery (%)
1	1 / 100	56.9	57.0	99.8
	1 / 200	27.6	28.5	96.8
	1 / 400	13.8	14.3	96.5
	1 / 800	6.9	7.1	97.2
2	1 / 100	98.6	100.0	98.6
	1 / 200	46.8	50.0	93.6
	1 / 400	24.9	25.0	99.6
	1 / 800	12.0	12.5	96.0

	Product Ref.	10215
	Product Desc.	β 2 GP1 IgGAM
	Manual Rev. No.	006: 2024-02-15

10.4 Precision

To determine the precision of the assay, the variability (intra and inter-assay) was assessed by examining its reproducibility on three serum samples selected to represent a range over the standard curve.

Intra-assay		
Sample No.	Mean (U/ml)	CV (%)
1	56.4	3.8
2	75.4	7.9
3	254.0	6.1

Inter-assay		
Sample No.	Mean (U/ml)	CV (%)
1	58.7	4.2
2	82.7	7.7
3	264.4	6.6

10.5 Calibration

Due to the lack of international reference calibration this assay is calibrated in arbitrary units (U/ml).

11 Literature

Schousboe, I. (1985): β 2-glycoprotein I: a plasma inhibitor of the contact activation of the intrinsic blood coagulation pathway. Blood, 66: 1086-1091.




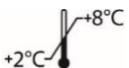

Nimpf, J, Bevers, EM, Bomans, PH, Till, U, Wurm, H, Kostner, GM, Zwaal, RF et al. (1986): Prothrombinase activity of human platelets is inhibited by β 2-glycoprotein I. Biochim Biophys Acta; 884: 142-149.

Nimpf, J, Wurm, H, Kostner, GM. (1987): β 2-glycoprotein I (apo-H) inhibits the release reaction of human platelets during ADP-induced aggregation. Arteriosclerosis; 6: 109-114.

Harris, E.N., Gharavi, A.E., Boey, M.L., et al. (1983): Anticardiolipin antibodies: Detection by radioimmunoassay and association with thrombosis in systemic lupus erythematosus. Lancet Nov 26, 1211-1214.

Galli, M, Comfurius, P, Maassen, C, et al. (1990): Anticardiolipin antibodies directed not to cardiolipin but to a plasma protein cofactor. Lancet 335: 1544-1547.

Wöhrle R, Matthias T, von Landenberg P, Oppermann M, Helmke K, Förger F (2000): Clinical relevance of antibodies against different phospholipids. Journal of Autoimmunity 15, A60.

IVD	- Diagnosi in vitro - Pour diagnostic in vitro - In Vitro Diagnostikum - Para uso Diagnóstico in vitro	- For in vitro diagnostic use - Para uso diagnóstico in vitro - In Vitro Διαγνωστικό μέσο
REF	* Numero d'ordine * Référence Catalogue * Bestellnummer * Número de catálogo	* Catalogue number * Numéro de catálogo * Αριθμός παραγγελίας
LOT	* Descrizione lotto * Lot * Chargen Bezeichnung * Lote	* Lot * Lote * Χαρακτηρισμός παρτίδας
UDI	* Identificatore univoco del dispositivo * Identifiant unique de l'appareil * eindeutige Produktidentifizierung * Identificador único do dispositivo	* Unique Device Identifier * Identificador único del dispositivo * Μοναδικό αναγνωριστικό συσκευής
CE	* Conformità europea * Déclaration CE de Conformité * Europäische Konformität * Declaração CE de Conformidade	* EC Declaration of Conformity * Declaración CE de Conformidad * Ευρωπαϊκή συμφωνία
	* 96 determinazioni * 96 tests * 96 Bestimmungen * 96 Testes	* 96 tests * 96 pruebas * 96 προσδιορισμοί
	* Rispettare le istruzioni per l'uso * Voir les instructions d'utilisation * Gebrauchsanweisung beachten * Ver as instruções de uso	* See instructions for use * Ver las instrucciones de uso * Λάβετε υπόψη τις οδηγίες χρήσης
	* Da utilizzarsi entro * Utilise avant le * Verwendbar bis * Utilizar antes de	* Use by * Utilizar antes de * Χρήση μέχρι
	* Conservare a 2-8°C * Conserver à 2-8°C * Lagerung bei 2-8°C * Conservar entre 2-8°C	* Store at 2-8°C (35.6-46.4°F) * Conservar a 2-8°C * Φυλάσσεται στους 2-8°C
	* Prodotto da * Fabriqué par * Hergestellt von * Fabricado por	* Manufactured by * Fabricado por * Κατασκευάζεται από
CON +	* Controllo positivo * Contrôle Positif * Positiv Kontrolle * Controllo positivo	* Positive Control * Control Positivo * Θετικός ορός ελέγχου
CON -	* Controllo negativo * Contrôle Négatif * Negativ Kontrolle * Controllo negativo	* Negative Control * Control Negativo * Αρνητικός ορός ελέγχου
CAL	* Calibratore * Etalon * Kalibrator * Calibrador	* Calibrator * Calibrador * Αντιδραστήριο βαθμονόμησης
CONJ	* Coniugato * Conjugué * Konjugat * Conjugado	* Conjugate * Conjugado * Σύζευγμα
MP	* Micropiastro rivestita * Microplaque sensibilisée * Beschichtete Mikrotiterplatte * Microplaca revestida	* Coated microtiter plate * Microplaca sensibilizada * Επικαλυμμένη μικροτράκα
WASHB 50x	* Tampone di lavaggio * Tampon de Lavage * Waschpuffer * Solução de lavagem	* Wash buffer * Solución de lavado * Ρυθμιστικό διάλυμα πλύσης
SUB	* Tampone substrato * Substrat * Substratpuffer * Substrato	* Substrate buffer * Tampón sustrato * Ρυθμιστικό διάλυμα υποστρώματος
STOP	* Reagente bloccante * Solution d'Arrêt * Stopreagenz * Solução de paragem	* Stop solution * Solución de parada * Αντιδραστήριο διακοπής αντίδρασης
SB 5x	* Tampone campione * Tampon Echantillons * Probenpuffer * Diluente de amostra	* Sample buffer * Tampón Muestras * Ρυθμιστικό διάλυμα δειγμάτων